# FYSIC<sup>®</sup>


# FW-150

Personenweegschaal

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Personenwaage

Bathroom scale

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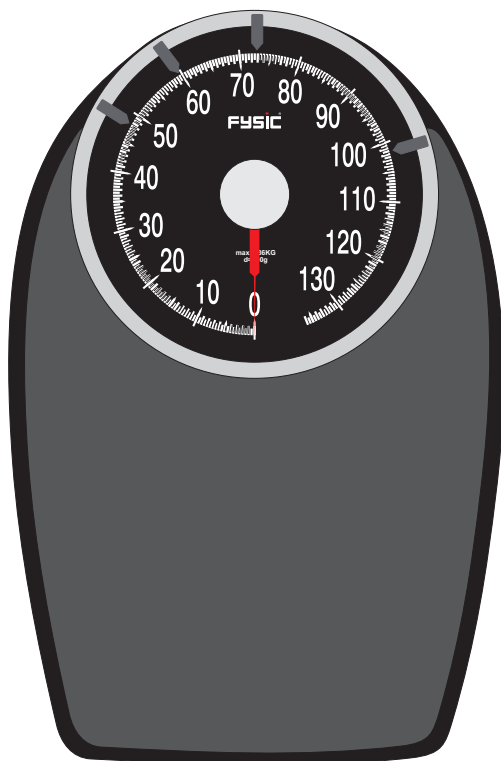Mode d'emploi

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User manual

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FR: Mode d'emploi: page 6

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
GB: User manual: page 10

## **TECHNICAL SPECIFICATIONS**

Max. weight:	136Kg.
Resolution:	500g.
Deviation 0-65Kg.:	$\pm 1,2$ Kg.
Deviation >65Kg.:	$\pm 2$ Kg.
Scale marker:	4 x

Geheugen marker  
marqueur de la mémoire  
Merker  
Scale marker

Gewicht  
Poids  
Gewicht  
Weight

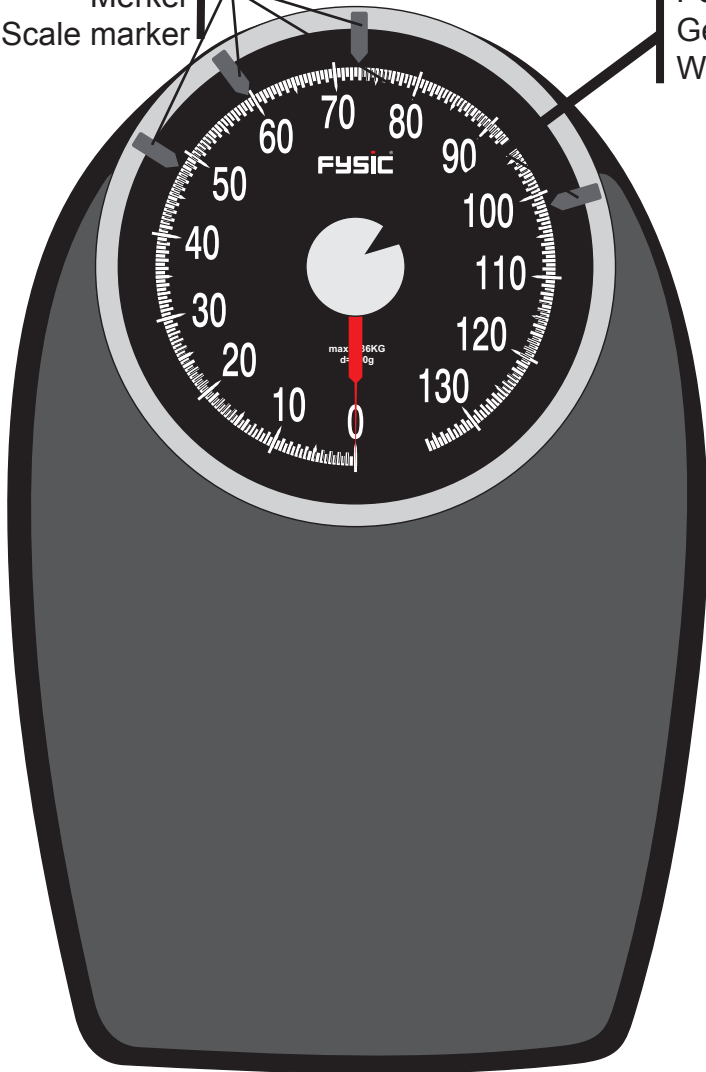0-Instelling  
0-Paramètre  
0-Einstellung  
0-Setting

## 1. BEOOGD GEBRUIK

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## 1. INTENDED USE

The Fysic FW-150 is a mechanical bathroom scale, suitable for weighing the weight of a person in a private household situation. The FW-150 is NOT suitable for hospital use or for industrial or commercial purposes.

## 2. IMPORTANT SAFETY INSTRUCTIONS

- Read the user manual carefully and completely, and make sure you understand all the instructions. It is also advised to keep this booklet as a reference.
- If you resale the scale, then ensure the booklet is with the device.
- The scale is not a toy. Please ensure that it is always kept away from children.
- It is advised to supervise the use of this scale together with persons with reduced physical or mental capabilities.
- Never immerse the scale in water.
- This scale is only for domestic use.
- Make sure the scale is placed horizontally on a firm, stable and level surface.
- Never place the scale on a flame or a hot surface .
- Ensure the floor under and around the scale is dry so as to avoid slipping.
- This scale does not contain any spare parts. Do not open the scale. Service repairs should only be performed by qualified technical personnel. Improper repairs could lead to a dangerous situation for the user.
- The FW-150 can only weigh up to 136Kg and NOT more.

### 3. USE

#### Preparation:

Make sure the scale is placed horizontally on a firm, stable and level surface.

Ensure that the floor around the scale is dry.

Take off any clothing in order to get a more accurate weight.

#### 0 setting:

Before anyone stands on the scale, use the scroll wheel located at the front to point the pointer exactly to 0Kg.

#### Weighing:

Step with both feet onto the middle of the scale. (Do not step too near the edge of the scale to prevent tipping). Read your weight shown by the pointer, and slide the memory marker on the measured value. 4 markers for 4 people.

### 4. CLEANING & MAINTENANCE

- Clean the outside of the scale with a damp cloth. Do not use any harsh or abrasive chemical cleaning agents.
- Apart from cleaning, the FW-150 does not require any further maintenance.
- The FW-150 does not contain any spare parts.


### 5. DISPOSING OF THE PRODUCT (ENVIRONMENTALLY)

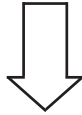When this product reaches the end of its lifecycle, do not throw this product into the normal household waste, rather take it to a collection point for the recycling of electrical and electronic equipment.

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